

THE TECHNOLOGY YOU WANT.
THE FLEXIBILITY YOU NEED.



The advanced technology of GentleYAG is more accessible than ever thanks to Pathways—your path to unprecedented Choice, Flexibility and Value. The **choice** of a complete range of innovative technologies for optimal treatment. The **flexibility** to meet your specific needs today, then easily upgrade as you expand your practice. The **value** that comes from cost-effectively growing your business—step-by-step, with the leader in medical and aesthetic lasers.

To learn more about Pathways, visit Candela at www.candelalaser.com/pathways.asp or call 1-800-821-2013.



Photographs courtesy of Mark B. Taylor, M.D., Laurie M. Hobbs, M.S., Simon Ourian, M.D., E. Victor Ross, M.D., Luigi L. Polla, M.D., Stephen W. Eubanks, M.D.

*Availability varies by country.

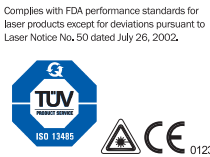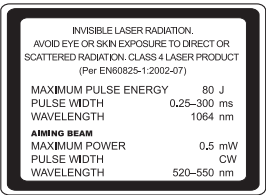

Specifications subject to change without notice.

GentleYAG and GentleYAG Limited Edition are registered trademarks of Candela Corporation. Dynamic Cooling Device and DCD are trademarks. ©2005 Candela Corporation. All rights reserved. Printed in the USA. 03/05 0920-11-0039 Rev. 05

SPECIFICATIONS



Laser Type	Nd:YAG laser
Wavelength	1064 nm
Repetition Rate	Up to 10 Hz
Pulse Duration	0.250–300 ms
Spot Sizes (diameter)	VR 1.5 mm, 3 mm SRT 6 mm, 8 mm, 10 mm HR 12 mm, 15 mm, 18 mm
Fluence Range	Up to 600 J/cm ²
Beam Delivery	lens-coupled optical fiber with handpiece
Pulse Control	Fingerswitch, Footswitch
Dimensions	35" H x 16" W x 28" D (89 cm x 41 cm x 71 cm)
Weight	210 lbs (95 kg)
Electrical	230 V 50/60 Hz, single phase, 16A

OPTIONAL DYNAMIC COOLING DEVICE

Integrated controls, cryogen container and handpiece with distance gauge	
Cryogen	HFC 134a
DCD Spray Duration	User adjustable range: 10-100 ms
DCD Delay Duration	User adjustable range: 3, 5, 10–150 ms
DCD Post Spray Duration	User adjustable range: 0–50 ms

Candela Corporation
530 Boston Post Road
Wayland, MA 01778 USA
Phone: (508) 358-7637
Fax: (508) 358-5569
Toll Free: (800) 821-2013
www.candelalaser.com



UNMATCHED Nd:YAG TECHNOLOGY



Build Your Own GentleYAG to Fit Your Practice and Budget

GentleYAG	✓	✓	✓	✓	✓	✓
GentleYAG HR	✓	✓		✓		✓
GentleYAG SRT	✓		✓			
GentleYAG VR		✓	✓			
GentleYAG HR/SRT	✓	✓		✓		✓
GentleYAG HR/VR	✓		✓		✓	
GentleYAG SRT/VR		✓	✓	✓		✓
GentleYAG LE	✓	✓	✓	✓	✓	✓

As your practice grows and diversifies, so should your laser. Candela's new GentleYAG is exclusively designed to grow with you. Multiple configurations—designed to treat a variety of conditions singularly or in combination—offer you the flexibility to tailor your practice to the needs of your patients.

Choose the Combination That's Right for You

GentleYAG® LEADING TECHNOLOGY FROM THE LEADER IN AESTHETIC LASERS



As a leader in aesthetic lasers since 1970, Candela is uniquely qualified to serve the changing needs of your growing laser practice. Our new GentleYAG comes in multiple configurations, features handpieces that are tailored to treatment applications, and patented cooling technology to offer you the most versatile Nd:YAG system ever made.



THE BEST FOR GROWING YOUR BUSINESS


Unlike other technologies that limit you to a one-size-fits-all approach, GentleYAG can treat virtually any cosmetic laser patient who walks through your practice door — efficiently, comfortably and profitably. How? GentleYAG's unique combination of a 1064 nm wavelength and Candela's patented Dynamic Cooling Device™ (DCD™) technology allow for greater protection of melanin-rich Fitzpatrick V and VI skin types, and for treating tanned skin. This enables you to offer comfortable, year-round procedures for the widest possible range of patient needs, including treatment of:

- Unwanted hair, for all skin types
- Skin Tightening
- Wrinkles
- Leg veins and facial veins
- Pseudofolliculitis Barbae (PFB)
- Vascular Lesions

A DYNAMIC COOLING DEVICE SO ORIGINAL, IT'S PATENTED

Candela's patented Dynamic Cooling Device™ (DCD™) provides consistent, reproducible epidermal cooling and unparalleled comfort for your patients. DCD prepares skin for treatment by selectively spraying the upper layers with a cooling burst of cryogen in adjustable spray/delay durations, milliseconds before and after the laser pulse. DCD's unique design provides excellent visibility throughout the entire process. The most powerful benefit to you? Speed and Patient Safety. DCD is operator technique independent—it will provide the same epidermal protection no matter how fast the laser is fired. The more patients you see and treat optimally, the more profitable for your practice..

with
DCD

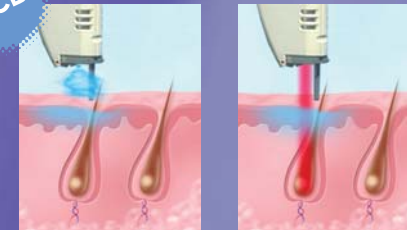

SUPPORT AFTER THE SALE

A GentleYAG purchase includes access to MyCandela.com, an exclusive marketing resource for your new laser. Services Include Website Builder, Laser marketing tool Kits. Learning center and additional clinical information on the capabilities of your Candela laser.*