

MULTIWell

THE SLIMMING CONCEPT FROM FRANCE



Ultra-sonic waves emission device **EXCEL CELL SONIC**

The low-frequency ultrasound with their lypolysis action penetrate the tissues and hammer the fat cells. Acting as micronisators, they break down fat particles. They also, at the same time, break up fibrinolytic partitioning and have a vasodilatory effect. The release and evacuation of fat is thus facilitated.

EXCEL CELL SONIC




- 4 independent generators of ultrasounds .
- Frequency : 30 kHz
varying between 28 and 32 kHz
- Power : 2W/cm² for each generator.
- 5 programmes of work.
- Cells of 1 to 3 transducers.
- Plane surface, polished and anodised.
- Water resistant cell and easy to clean.
- Protection against over-heating.



MULTIWELL
ZI DU GRAND PONT
13640 LA ROQUE D'ANTHERON
Tl. : 04.42.28.48.20
Fax. : 04.42.28.50.81

MULTIWell

THE SLIMMING CONCEPT FROM FRANCE

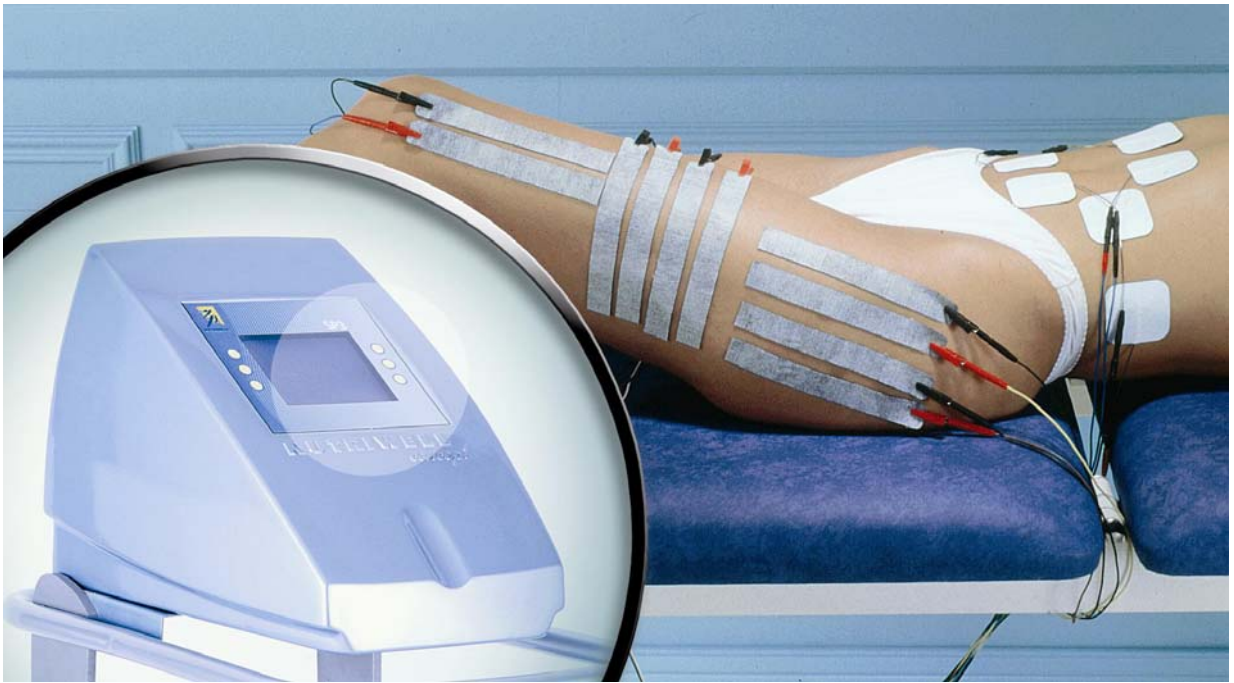

Electrostimulation device **EXCEL SP 3**

The SP 3 is an extremely high performance stimulator. Its microprocessors and 12 outlets enable it to simultaneously manage different stimulation programmes. The SP 3 is equipped with specific programmes for each part of the body or each type of use. You just select the right programme and adjust the intensity level to that required for the treatment. **The SP 3 generators are equipped with securised outlets to protect the customer from any risk of electric shock or burning.**

EXCEL SP 3



Safe outputs : each output is safe, when an electrode is unplugged, There is no current

Electrical characteristics

Power supply : 220Volts / 50Hz (alternative ; sinusoidal)

Power consumption : 64W

Apparatus of type IIBF

Fuses : 2 fast fuses 240V / 0.3A 4/20mm

In accordance with directive 93/42/CEE



MULTIWELL
ZI DU GRAND PONT
13640 LA ROQUE D'ANTHERON
Tl. : 04.42.28.48.20
Fax. : 04.42.28.50.81

MULTIWell

THE SLIMMING CONCEPT FROM FRANCE



***DLS 10** Pressure-therapy device*

Pressoplastie promotes the transfer of fat molecules in the lymphatic system to the venous system and then on to the burners. Cell regeneration induced by DLS 10 is similar to that produced by manual drainage. DLS 10 is equipped with a microprocessor to enable the management of series of automatic programmes, extremely easy to use and very high performance.

***DLS 10** Pressure-therapy device*



General Characteristics

- Principle : Sequential pressure applied : electromagnetic value 20
- Number of outlets : 3x10 tube couplers with central locking.
- Sealing : Couplers may be sealed with a plug.
- Sequence time : 10 seconds, adjustable from 4 to 20 seconds in stages of 1 sec.
- Pressure : 0-150 mm Hg in stages of 1mmHg. Regulated by incorporated pressure sensor.
- Session time : 1 to 60 minutes in stages of 1 minute.
- Display : large backlit screen.
- Mains Supply : 220V – 50/60 Hz.



MULTIWELL

ZI DU GRAND PONT


13640 LA ROQUE D'ANTHERON

Tl. : 04.42.28.48.20

Fax. : 04.42.28.50.81

MULTIWell

THE SLIMMING CONCEPT FROM FRANCE

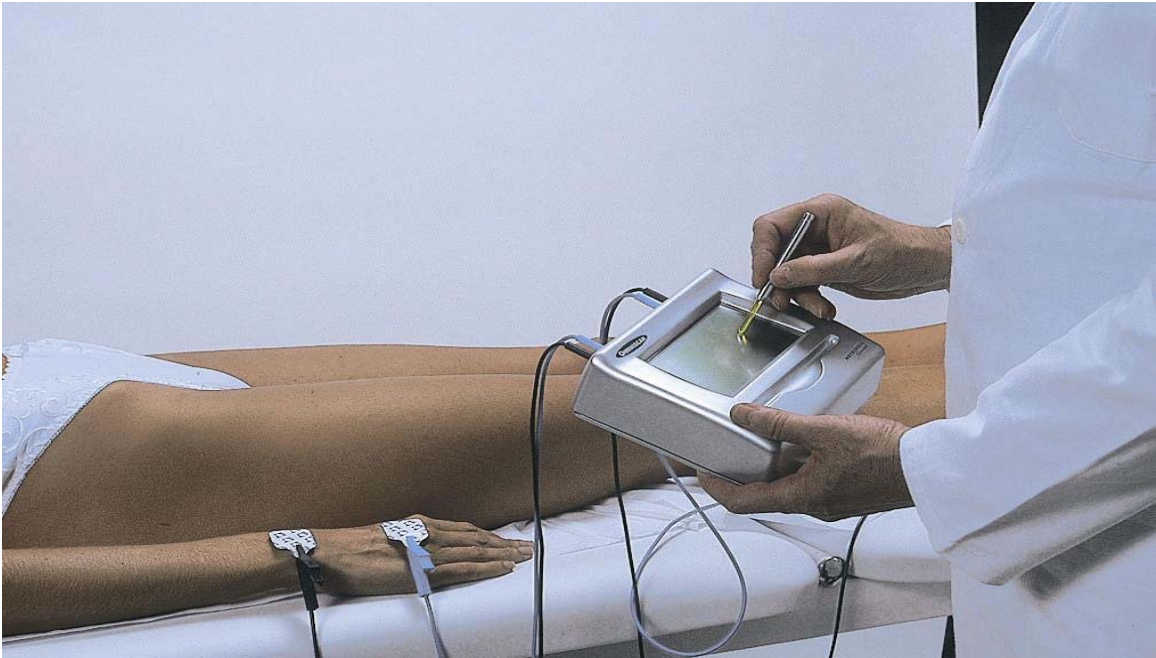

OHMNISCAN NG

This method is based on the conduction of an alternative 50 kHz current, of very low intensity – 400 μ A – through the entire body. Organic tissue acts as an ionic conductor, that is to say that in order to circulate in this type of conductor, the current uses ions in suspension in the body fluid. The greater the quantity of water, or the more free ions there are, the more easily the current circulates. Measurement thus depends on total fluid. Using impedance analysis – resistance of the body to the passage of current – the body's TOTAL FLUID is measured by the relationship illustrated opposite.


OHMNISCAN NG



General Characteristics

- Type of measurement: 4-wire differential.
- Measurement range : 100 – 1,500 ohms.
- Measurement precision : 2% at 500 Ω at 25°C throughout the life of the batteries
- Generator : sinusoidal 50 kHz – 0,4 μ A efficiency.

METHOD OF MEASUREMENT

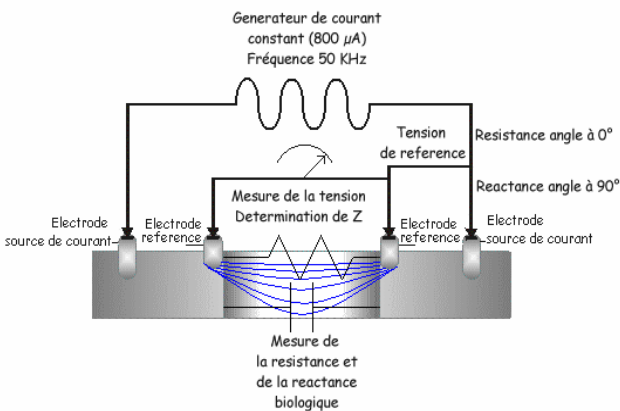

- Communication : infra-red link (irda 2400 bds).
- Supply: 4 CdNi 1,5 V batteries screen and battery charging with special medical transformer.
- Battery charge : approx. 300 hours.



MULTIWELL
ZI DU GRAND PONT
13640 LA ROQUE D'ANTHERON
Tel. : 04.42.28.48.20
Fax. : 04.42.28.50.81

MULTIWell

THE SLIMMING CONCEPT FROM FRANCE



Vacuum therapy device **Bioderm NG**

The Bioderm NG is a combination of two techniques : cupping and skin-rolling. The combination of these techniques as part of a holistic Slimming concept increases tenfold the results for any weight problem, Including the most stubborn cases.

Its transparent heads, enables skin folds to be formed through aspiration, which can then be moved over all the parts of the body. Cellulite can be removed at each point of formation and general circulation in the upper layers of the epidermis can be improved.

BIODERM NG



- Hyper-vascularisation
- A regulation of the Blood flow
- A break-down of the collagenous fibers
- A relaxing and anti-stress effect
- A draining effect

Specifications

- 50 to 800 mBar in 15 stages
- 3 body drainage
- 3 slimming programmes
- 4 face care programmes
- 3 scar treatment programmes
- Input voltage: 220 /240 v –50/60 Hz



MULTIWELL

ZI DU GRAND PONT

13640 LA ROQUE
D'ANTHERON

Tel. : 04.42.28.48.20

Fax. : 04.42.28.50.81